

## ELISAPro Validation: Small Ruminant Lentivirus (SRLV) Antibody Test Kit, cELISA

Each assay undergoes the following series of tests:


**Timing** Confirms testing is completed according to assay circular, timing based per sample with a well-to-well precision of < 15 seconds

**Result:**

	Well 1	Well 48	Well 88
Sample	1:00:10	1:00:05	1:00:05
Conjugate	0:30:08	0:30:07	0:30:06
Substrate	20:08	20:08	20:09

**ELISAPro vs Manual** Compares the data from ELISAPro to data obtained manually from two technicians for confirmation of concordance using an internal, well-characterized sample set

**Result:**

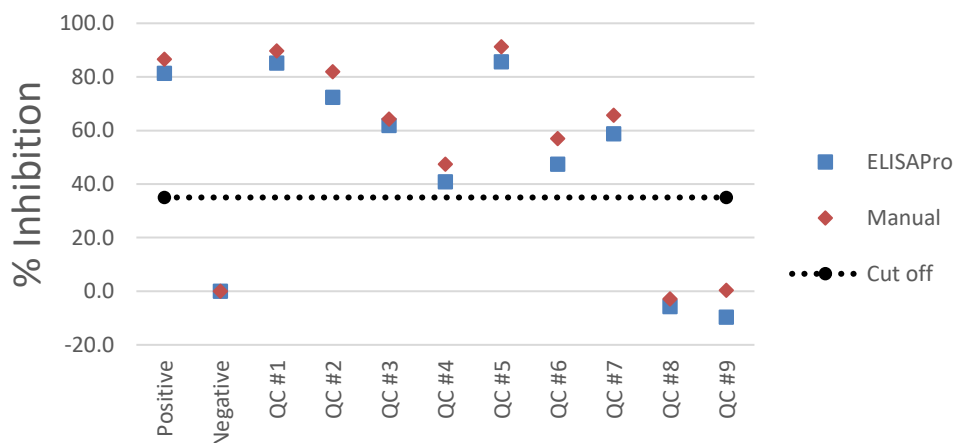

**Cross-contamination** Confirms there is no well to well cross-contamination due to carryover of samples

**Result:**

The assay was negative for cross-contamination (data not shown).

**Conclusion**

In accordance with OIE standards, VMRD has demonstrated that the ELISAPro adequately performs VMRD's SRLV Antibody Test Kit, cELISA and produces data that is equivalent to that obtained using traditional methods using a well-characterized panel of samples. Each sample and control met pre-determined criteria for both methods.