


CERTIFICATE OF ANALYSIS

## Canine Herpesvirus Type 1 (CHV-1)

FITC Conjugate

Catalog No.:	CJ-F-CHV-10ML
Volume:	10 ml
Lot:	P140407-003
Expiration:	03 January 2022
Agent:	Canine Herpesvirus Type 1 (CHV-1)
Strain:	NA

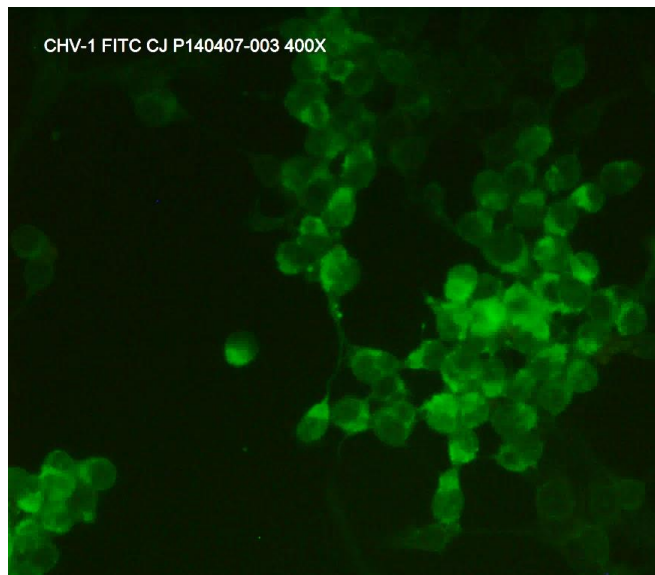

### Description:

CHV-1 polyclonal antiserum conjugated to fluorescein isothiocyanate (FITC). Canine origin. Ready to use. Liquid.

### Quality Control Method:

Direct FA using CHV-1 12 well slide (catalog no. SLD-IFA-CHV).

**Specific Reaction:** 2-3+ signal with trace background.

**Other Comments:** The raw material has also been screened by Indirect FA and does not react with canine adenovirus type 1 and 2 (CAV-1 and 2), canine coronavirus (CCV), canine distemper virus (CDV), canine parainfluenza virus type 2 (CPI-2), canine parvovirus (CPV), reovirus (REO), Rabies Recombinant Nucleoprotein (rNP), *Toxoplasma gondii*, vesicular stomatitis virus Indiana and New Jersey strains (VSV).

### Pattern Of Fluorescence:

Cytoplasmic and nuclear fluorescence in rounded cells with occasional small syncytia.

### Intended Use:

This reagent is suitable for CHV-1 virus identification in cell cultures and in animal tissues.

### Storage:

This conjugate is provided in liquid form and should be stored at 2-7°C. DO NOT FREEZE! It should also be stored in the original container and/or in the dark. If conjugate becomes cloudy it should be discarded. This conjugate contains 0.09% sodium azide as a preservative.

### References:

Ledbetter EC, et al., Experimental reactivation of latent canine herpesvirus-1 and induction of recurrent ocular disease in adult dogs. Vet Microbiol. 2009 Jul 2; 138(1-2):98-105. Epub 2009 Mar 13.

#### Recommended Staining Procedure for Indirect FA:

1. Warm slide to room temperature before removing from foil pouch.
2. Place diluted serum on the designated wells. Dilute serum in serum diluting buffer, pH 7.2 (catalog no. FASDB-100ML) however if high background due to anti-bovine IgG activity is present it may be advisable to use SSDB-100ML.
3. Incubate slide in humid chamber at 37°C for 30 minutes.
4. Using a wash bottle, gently rinse slide briefly in FA rinse buffer, pH 9.0 (catalog no. FARB-4X) and then soak for 10 minutes in FA rinse buffer, pH 9.0.
5. Drain slide and dry around wells by pressing blotter (included in pouch) to front surface. Place labeled anti-IgG or IgM on the wells.
6. Incubate as in step 3.
7. Rinse as in step 4.
8. Drain slide and dry back and edges with a paper towel. Do not allow stained surface to dry. Do not rinse with water.
9. Mount with mounting fluid [glycerol/FA rinse buffer, pH 9.0, (50/50)] (catalog no. FAMF-10ML) and view with good quality fluorescence microscope at 100X-250X. Confirmation may be made at 400X.

#### Recommended Staining Procedure for Direct FA:

1. Warm slide to room temperature before removing from foil pouch.
2. Place direct FA conjugate on the designated wells.
3. Incubate slide in humid chamber at 37°C for 30 minutes.
4. Using a wash bottle, gently rinse slide briefly in FA rinse buffer, pH 9.0 (catalog no. FARB-4X) and then soak for 10 minutes in FA rinse buffer, pH 9.0.
5. Drain slide and dry back and edges with a paper towel. Do not allow stained surface to dry. Do not rinse with water.
6. Mount with mounting fluid [glycerol/FA rinse buffer, pH 9.0, (50/50)] (catalog no. FAMF-10ML) and view with good quality fluorescence microscope at 100X-250X. Confirmation may be at 400X.

#### Serum Diluting Buffer (pH 7.2):\*

- Na<sub>2</sub>HPO<sub>4</sub>.....1.19 gm
- NaH<sub>2</sub>PO<sub>4</sub>.....0.22 gm
- NaCl.....8.55 gm
- BSA.....10.0 gm
- DI/dH<sub>2</sub>O.....Q.S. to 1 liter

\*This recipe makes 1 liter. If you need less, adjust recipe accordingly. Store at 2-7°C. Add 0.09% NaN<sub>3</sub> if diluted serum is not going to be used within one week.

#### 4X FA Rinse Buffer (pH 9.0):

- Na<sub>2</sub>CO<sub>3</sub>.....11.4 gm
- NaHCO<sub>3</sub>.....33.6 gm
- NaCl.....8.5 gm
- DI/dH<sub>2</sub>O.....Q.S. to 1 liter

Final pH should be 9.0-9.5. This is a 4X concentrate and should be diluted 1/4 with DI/distilled water for use as a working buffer. Keep in a tightly stoppered container at room temperature. MOUNTING FLUID is made by mixing glycerol and FA rinse buffer, pH 9.0, in equal proportions.