


Equine Infectious Anemia Virus (EIAV), ELISA Medium Positive

Reference Serum

Catalog No.:	RS-EIA-EM-0.5ML
Shelf Life:	2 years from date of qualification
Volume:	0.5 ml
Immunogen:	Equine Infectious Anemia Virus (EIAV)
Preservative:	10 ppm ProClin 300
Species of Origin:	Equine

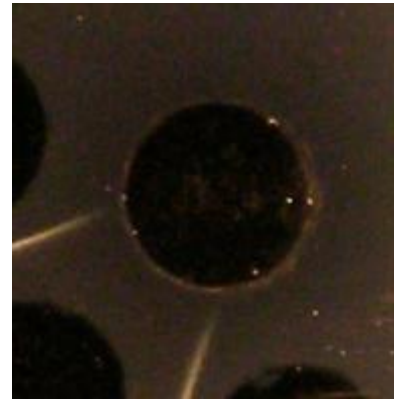

Description:

Equine serum containing a level of antibody to Equine Infectious Anemia Virus (EIAV) that gives a medium positive reaction in the VMRD ELISA and a weak positive in the VMRD AGID.

Intended Use:

This product is intended as a reference sample for quality assurance of EIA ELISA tests.

Interpretation of Results:

The sample's expected OD range in the VMRD EIA ELISA is > 0.25 and < 0.75 . It is also usually weakly positive by AGID.

Storage:

This reference serum is provided in liquid form and should be stored at 2-7°C. If serum becomes cloudy, it should be discarded.

References:

NA

Technical Data Sheet Version:

Version 1

FOR IN VITRO LABORATORY USE ONLY.

H:\Master Documents\TDS\Reference Serum\RS-EIA-EM-0.5ML TDS V1 Active MKT.docx

25 September 2014

