


**CERTIFICATE OF ANALYSIS**

*Anaplasma phagocytophilum*

FA Substrate Slide

Catalog No.:	SLD-IFA-AP
Size:	12 Well
Well Capacity:	10 µl
Lot:	P140205-003
Expiration:	10 February 2018
Agent:	<i>Anaplasma phagocytophilum</i>
Strain:	Martin
Cell Culture Substrate:	HL-60 cells

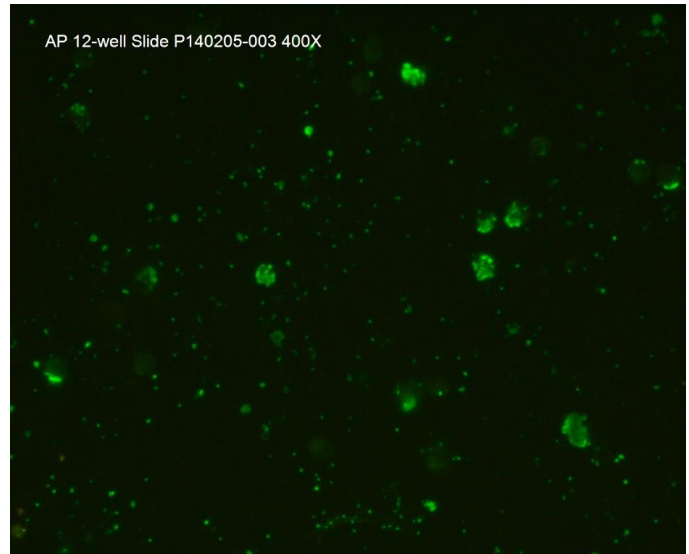

**Description:**

Slides contain fixed *A. phagocytophila* in HL-60 cells. Slides are unstained and sealed in moisture-free foil pouches.

**Quality Control Method:**

Indirect FA using *Anaplasma phagocytophilum* positive control (catalog no. PC-IFA-AP), negative control (catalog no. NC-IFA-AP), anti-equine IgG FITC conjugate (catalog no. CJ-F-EQUG-1ML or 10ML).

**Specific Reaction:** 3-4+ fluorescence with the positive control, no background. The negative control was negative with no background. There are 10-40 positive cells per high-power field.

**Other Comments:** NA

**Pattern Of Fluorescence:**

Brightly fluorescent cytoplasmic inclusion bodies (morulae), clustered elementary bodies (initial bodies), and free elementary bodies.

**Interpretation Of Results:**

Titers of 1/50 (IgG) and greater are considered antibody positive. Sera positive at the 1/50 screening dilution should be rerun to determine their endpoint titer for comparison with earlier or later specimens from the same animal.

**Intended Use:**

For detection of antibody to *A. phagocytophilum* by indirect FA technique. May be used to differentiate antibody class (IgG or IgM) with suitable quality fluoresceinated second antibody conjugate.

**Storage:**

Store sealed in foil pouch at -20°C. Avoid self-defrosting freezers.

**References:** NA

Recommended Staining Procedure for Indirect FA:

1. Warm slide to room temperature before removing from foil pouch.
2. Place diluted serum on the designated wells. Dilute serum in serum diluting buffer, pH 7.2 (catalog no. FASDB-100ML) however if high background due to anti-bovine IgG activity is present it may be advisable to use SSDB-100ML.
3. Incubate slide in humid chamber at 37°C for 30 minutes.
4. Using a wash bottle, gently rinse slide briefly in FA rinse buffer, pH 9.0 (catalog no. FARB-4X) and then soak for 10 minutes in FA rinse buffer, pH 9.0.
5. Drain slide and dry around wells by pressing blotter (included in pouch) to front surface. Place labeled anti-IgG or IgM on the wells.
6. Incubate as in step 3.
7. Rinse as in step 4.
8. Drain slide and dry back and edges with a paper towel. Do not allow stained surface to dry. Do not rinse with water.
9. Mount with mounting fluid [glycerol/FA rinse buffer, pH 9.0, (50/50)] (catalog no. FAMF-10ML) and view with good quality fluorescence microscope at 100X-250X. Confirmation may be made at 400X.

Recommended Staining Procedure for Direct FA:

1. Warm slide to room temperature before removing from foil pouch.
2. Place direct FA conjugate on the designated wells.
3. Incubate slide in humid chamber at 37°C for 30 minutes.
4. Using a wash bottle, gently rinse slide briefly in FA rinse buffer, pH 9.0 (catalog no. FARB-4X) and then soak for 10 minutes in FA rinse buffer, pH 9.0.
5. Drain slide and dry back and edges with a paper towel. Do not allow stained surface to dry. Do not rinse with water.
6. Mount with mounting fluid [glycerol/FA rinse buffer, pH 9.0, (50/50)] (catalog no. FAMF-10ML) and view with good quality fluorescence microscope at 100X-250X. Confirmation may be at 400X.

Serum Diluting Buffer (pH 7.2):\*

- Na<sub>2</sub>HPO<sub>4</sub>.....1.19 gm
- NaH<sub>2</sub>PO<sub>4</sub>.....0.22 gm
- NaCl.....8.55 gm
- BSA.....10.0 gm
- DI/dH<sub>2</sub>O.....Q.S. to 1 liter

\*This recipe makes 1 liter. If you need less, adjust recipe accordingly. Store at 2-7°C. Add 0.09% NaN<sub>3</sub> if diluted serum is not going to be used within one week.

4X FA Rinse Buffer (pH 9.0):

- Na<sub>2</sub>CO<sub>3</sub>.....11.4 gm
- NaHCO<sub>3</sub>.....33.6 gm
- NaCl.....8.5 gm
- DI/dH<sub>2</sub>O.....Q.S. to 1 liter

Final pH should be 9.0-9.5. This is a 4X concentrate and should be diluted 1/4 with DI/distilled water for use as a working buffer. Keep in a tightly stoppered container at room temperature. MOUNTING FLUID is made by mixing glycerol and FA rinse buffer, pH 9.0, in equal proportions.