



VMRD, Inc.

PO Box 502, Pullman, WA 99163 USA

Telephone: + 1 (509) 334-5815

Fax: + 1 (509) 332-5356

E-mail: vtech@vmrd.com

Web site: <http://www.vmrd.com>

Certificate of Analysis


Primary Anti-Feline Leukemia Virus (FeLV) Antiserum for IFA

CATALOG NO.: 210-46-FELV1

VOLUME: 10 ml

LOT: V02120-102504

EXPIRATION: 25 October 2004

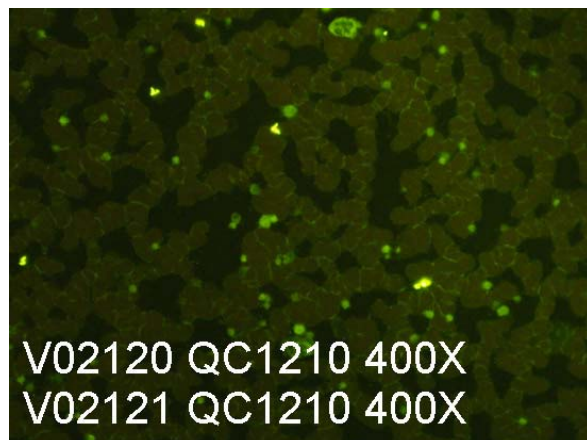

DESCRIPTION: Anti-FeLV polyclonal antiserum. Goat origin. Useful for detection of Feline Leukemia Virus by IFA (procedure follows). Anti-FeLV antiserum specifically binds FeLV antigens in peripheral blood leukocytes, necropsy tissues, or tissue culture.

QUALITY CONTROL METHOD: IFA on FeLV-infected blood smears using VMRD, Inc. FeLV positive blood smear slides (catalog no. 214-P-FELV), FeLV negative blood smear slides (catalog no. 214-N-FELV), FITC-labeled anti-IgG FITC Conjugate (catalog no. 210-47-FELV2).

Specific Reaction: 2-4+ on positive blood smear; negative on negative blood smear.

Other Reactions or Comments: Trace background (some artifacts).

STORAGE: Store at 4-8°C. If the antiserum becomes cloudy, it should be discarded. This antiserum contains 0.09% sodium azide as a preservative.

FOR RESEARCH AND INVESTIGATIONAL USE ONLY.

WARRANTY: VMRD, Inc. warrants that this product is as described in the quantity and contents stated on the label at the time of delivery to the customer. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE MADE BEYOND THE LABEL DESCRIPTION, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. Remedy is limited to replacement of the product or refund of the purchase price. VMRD, Inc. is not liable for property damage, personal injury, or economic loss caused by the product. The information listed in this information sheet is provided for reference only, and should not be substituted for the user's own incoming material quality control.

Detection of Feline Leukemia Virus by Indirect Fluorescent Antibody

I. PERIPHERAL BLOOD LEUKOCYTES

1. Buffy coat smears are preferred, but thin, feathered-edge blood smears (as made for a WBC differential) may be used for this test. They should be air dried, but not over one week old. Unstained smears should be stored at room temperature in a dark, dry place and should NOT be refrigerated due to the fact that any moisture condensing on the surface of the slide will ruin it for staining. Unstained smears may be preserved long term by storing in acetone-methanol (75/25%) at Revco temperatures (-76°C to -90°C). Fresh blood in anticoagulant should be mixed very thoroughly before smearing on clean, new slides. Fresh slides should be dried at least 30 minutes at room temperature or 15 minutes at 37°C before fixing.
2. Slides should be fixed for 20 minutes at room temperature in a mixture of 75% acetone and 25% methanol. This fixative may be mixed up in bulk in advance and stored at room temperature in a tightly capped vessel. Allow the slides to air dry after fixing and then circle (10-15 mm) the feathered edge with a Mark-Tex pen, liquid embroidery, grease pencil, or other suitable means.
3. Place 50-100 µl (1-2 drops) of Anti-FeLV Antiserum (VMRD catalog no. 210-46-FELV1), white label, primary antibody) in the circle and gently spread to cover. Incubate slides 30 minutes at 37°C in a humid chamber.
4. At the end of the primary incubation period, rinse the slide gently for 5-10 seconds in FA Rinse Buffer (VMRD catalog no. 210-90-RB) and then soak in FA rinse buffer for 10 minutes. Drain the slide and dry around the circle with a Kleenex or paper towel. Do not allow the smear inside the circle to dry.
5. Place 1-2 drops of Anti-IgG FITC Conjugate (VMRD catalog no. 210-47-FELV2, yellow label, secondary fluoresceinated antibody) in the circle and incubate 30 minutes at 37°C in a humid chamber.
6. Rinse as before for 10 minutes and drain. Dry the *back* of the slide (do not allow the front to dry).
7. Place one drop of mounting medium (50% FA rinse buffer/50% glycerol; VMRD catalog no. 210-92-MF) on the circle, mount with a coverslip and read on a good quality fluorescence microscope. Scan at 100-250X and confirm at 400X.
8. Stained RBC's will appear red to orange/brown. FeLV negative leukocytes will not be visible except for eosinophils, which will be white to pale yellow. Positive cells will appear bright apple green: neutrophils, eosinophils, lymphocytes, platelets and/or large lymphoblasts. The cells may be concentrated at the edge of the feathered area and large clumps of positive platelets or lymphocytes may be seen.

II. POST MORTEM TISSUES:

1. Impression smears of the cut surface of spleen, thymus, or lymph node may be made and processed as for blood smears. Frozen sections of any tissue may also be used and stained as above.
2. Bone marrow aspirates may be stained either from live cats or from post-mortem sources. Impressions or sections of biopsied tumors may be stained as an aid in the identification of FeLV-induced lymphosarcomas.

III. TISSUE CULTURES:

1. FeLV infection in tissue cultures may also be confirmed by this method. The protocol in Section I should be used, with the substitution of pure (100%) acetone fixation for 10 minutes instead of 75% acetone and 25% methanol as in Step 2 of Section I.
2. This protocol may be used in performing serum neutralization tests on cat sera by staining the culture 5-7 days after planting.

4X FA RINSE BUFFER (pH 9.0):

- Na ₂ CO ₃	11.4 gm	Final pH should be 9.0-9.5. This is a 4X concentrate and should be diluted 1:4 with DI/distilled water for use as a working buffer. Keep in a tightly stoppered container at room temperature.
- NaHCO ₃	33.6 gm	
- NaCl	8.5 gm	
- DI/dH ₂ O	Q.S. to 1 liter	